

HCL0032_HCL0033_HCL0034_PI_V001

Progesterone Receptor Analyte Control^{DR}

Cell line controls for immunohistochemistry
and in situ hybridization.

Research Use Only

PRODUCT AVAILABILITY

Product Code	Product Description
HCL0032	X2 Cut slides
HCL0033	X5 Cut slides
HCL0034	X1 Cell microarray block

APPLICATION

This product is suitable for use in immunohistochemistry.

MATERIALS

Four formalin fixed paraffin embedded cell lines with a dynamic range (DR) of expression for Progesterone Receptor (PR).

Cell line A: Ductal Carcinoma


Cell line B: Breast Adenocarcinoma

Cell line C: Ductal Carcinoma


Cell line D: Ductal Carcinoma

Cells are fixed in 10% neutral buffered formalin and paraffin wax embedded.

Sections are cut at 4µm, mounted on positively charged slides and baked overnight at 37°C.

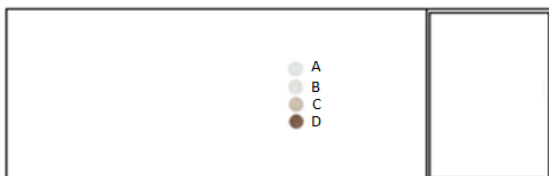

Cell microarrays (CMA) contain cores that are 2mm in diameter and 3-3.5mm in length. It is possible to obtain over 300 sections depending on thickness.

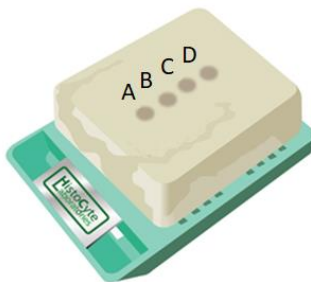

Expression Profile

Cell Line	IHC for Progesterone Receptor
A	Negative
B	Low/Intermediate
C	Intermediate/High
D	High

Storage and Handling

Store at 2-8°C. Do not freeze (for expiration date please see the product label)

WARNINGS AND PRECAUTIONS

1. The product is intended for research use only.
2. It is the responsibility of the end user to determine suitability with their reagents and procedures within their laboratory.
3. Do not use after expiration date printed on product labels. The user must validate any storage conditions other than those specified in the package insert.

TROUBLE SHOOTING

For further help please feel free to contact HistoCyte Laboratories Ltd at info@histocyte.com or call +44 (0)191 603 1007.



For updates and additional product information please visit: www.HistoCyte.com